



Read readme.md to run the server then go to <http://localhost:3000>

1. Upload Video Files

Click the Upload Video button in the left panel and select your video file(s) to upload.

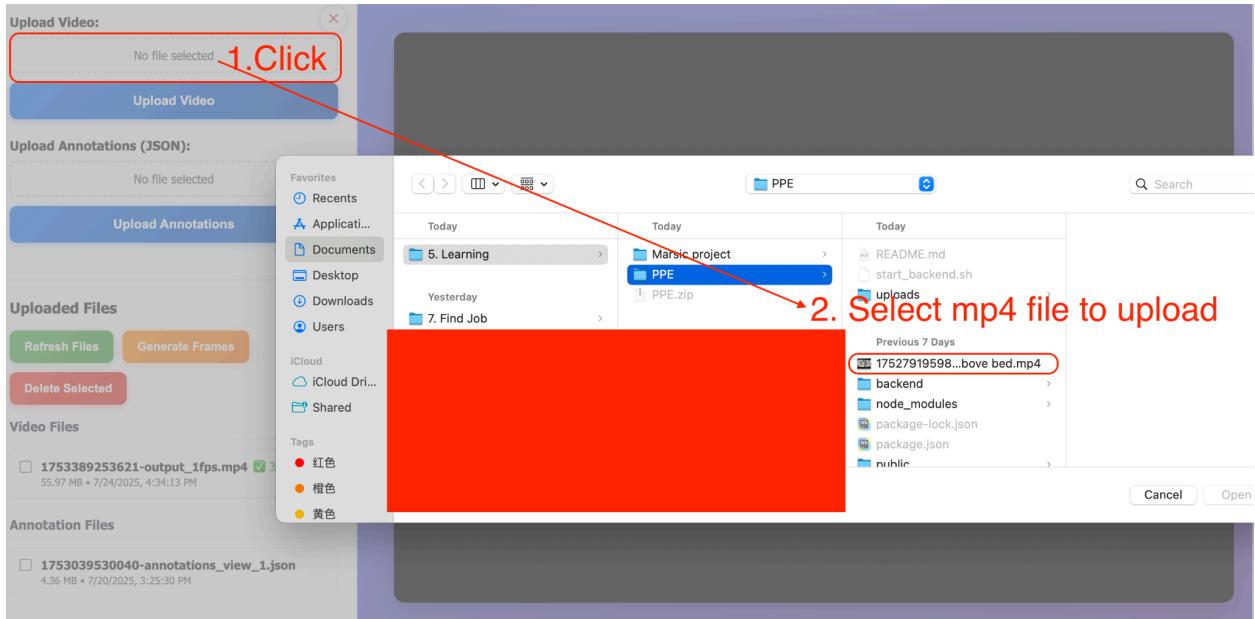


Wait for the upload and frame extraction to complete. Progress will be shown in the interface.

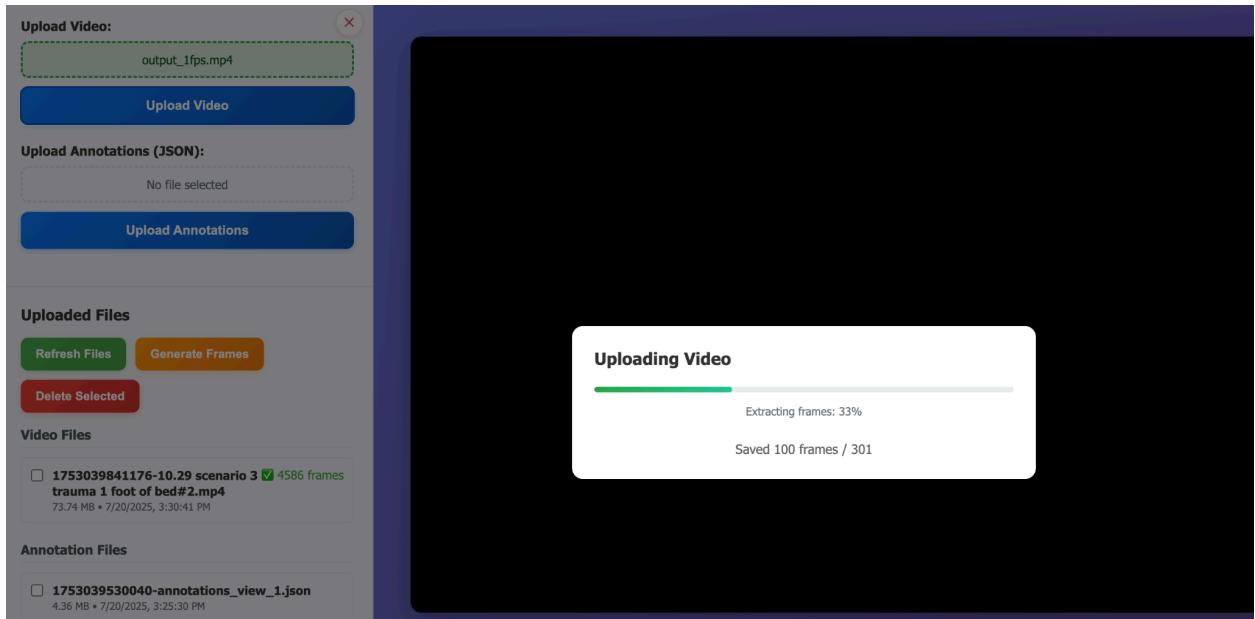

2. Upload Annotation Files

Click the Upload Annotation button and select the corresponding annotation JSON file for your video.

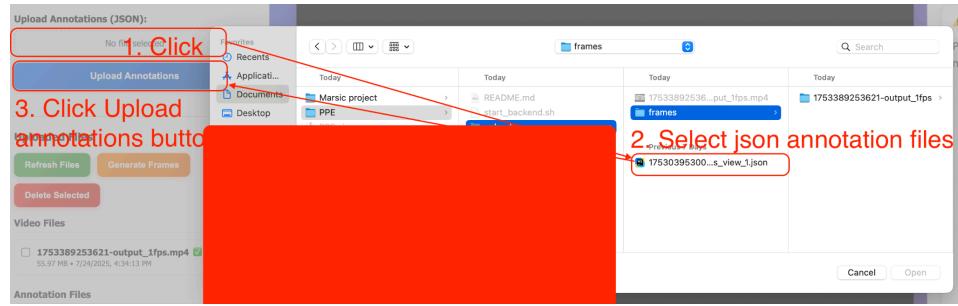

3. View and Annotate

Select a video and json annotation file from the file list.

1. Click an mp4 file
2. Click a json annotation file

Use the video controls to play, pause, or navigate frame by frame. Bounding boxes and PPE compliance information will be displayed on the video canvas.


1. Progress bar and play/pause button to control the video

Compliance Rate
16.7% (1/6)


PPE Uncompliance Details

- Member ID: 1
 - Gown Improper
 - Mask Absent
- Member ID: 3
 - Gown Absent
 - Gloves Absent

Click inside a bounding box to view detailed PPE non-compliance information in the right panel.

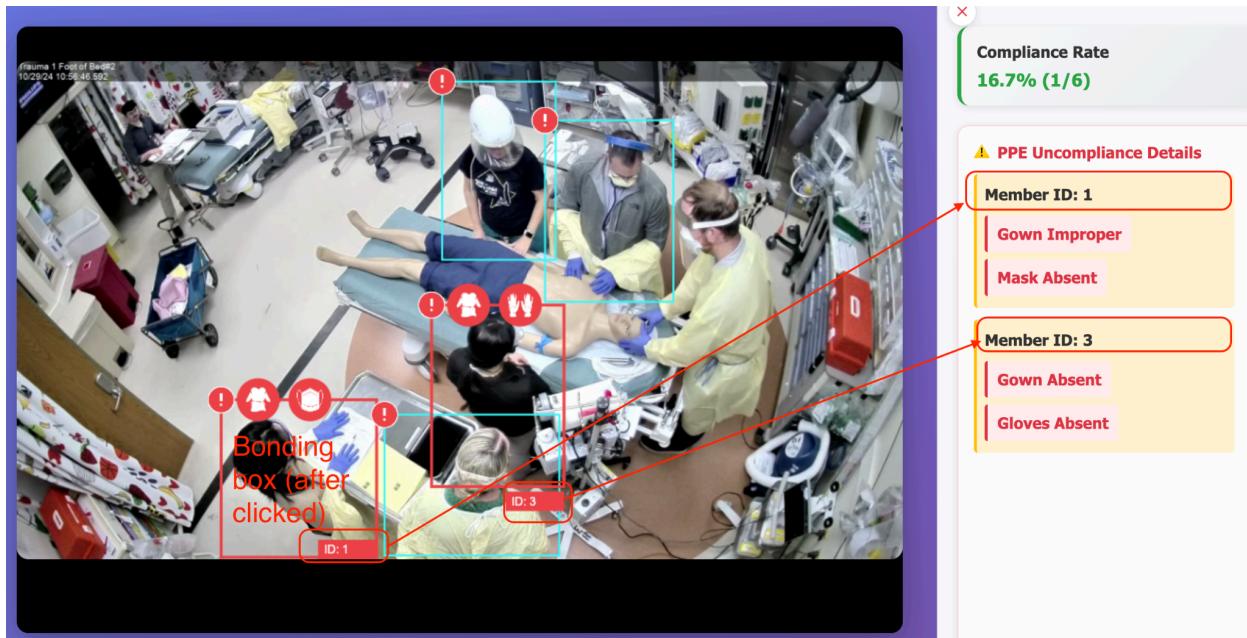

Click the **Show all** button to see all the bonding box in the current frame.

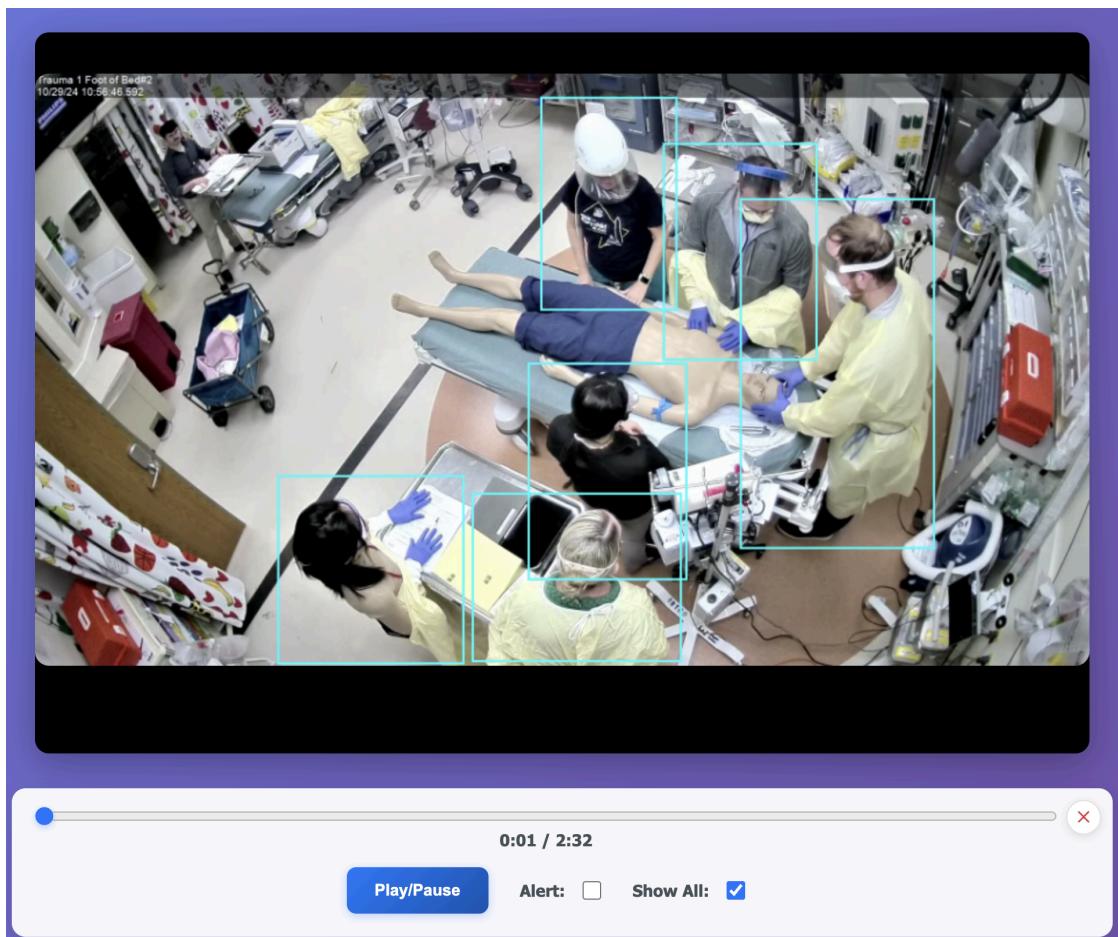


4. Review PPE Compliance

The right panel displays PPE non-compliance details for all selected individuals. You can select multiple bounding boxes to compare PPE status. The compliance rate and detailed missing PPE items are shown and updated in real time.

5. Additional Features

Hide or show the video controls bar using the triangle button at the bottom center.

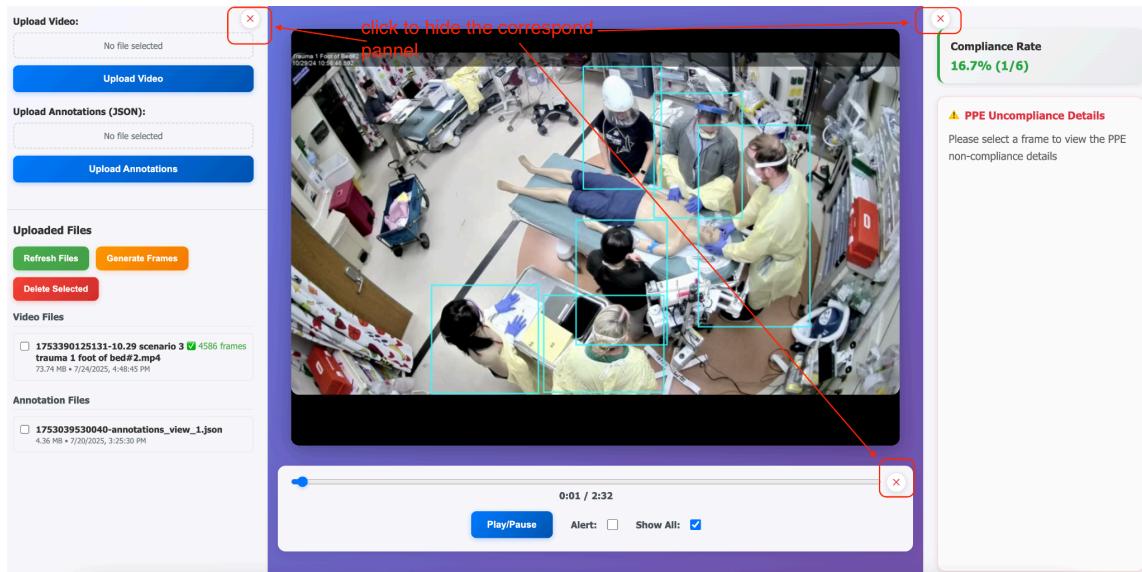

Use the file management buttons to refresh, generate frames, or delete files as needed.

If frames are not automatically extracted, select the video and click Generate Frames. The system will process the video and extract frames for annotation.

